



GelPOINT[®]

V-Path

Transvaginal Access Platform

GelPOINT V-Path

Transvaginal Access Platform

vNOTES-Enabling Technology

Provides vaginal access for minimally invasive instruments, enabling surgeons to perform vNOTES (vaginal natural orifice transluminal endoscopic surgery).





GelSeal® Cap

Provides a flexible fulcrum for excellent triangulation of standard laparoscopic instrumentation.

Includes two valves to facilitate smoke evacuation and accommodate standard laparoscopic insufflation.



Alexis® O Retractor

Enables transvaginal access by providing circumferential retraction.

Facilitates the extraction of large specimens.

Comes in three sizes to meet procedural needs.



Introducer for the Alexis O Retractor

Facilitates placement of the retractor into the colpotomy.


Comes in two designs to meet procedural needs.



Self-Retaining Sleeves

Accommodate 5–12mm standard laparoscopic instrumentation.

Allow for ample freedom of movement and reduce instrument collision as a result of their low-profile design.

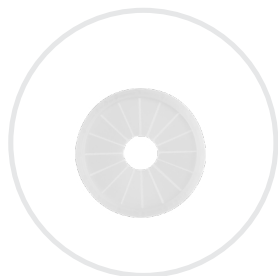

Optional Instrument Shield

Guides instruments through the incision.

GelPOINT V-Path Transvaginal Access Platform



| Reorder No. | Alexis O Retractor Inner Ring Size | Qty/Box |
|-------------|------------------------------------|---------|
| C2A11 | 11cm | 1 |

Components

| |
|---|
| (1) GelSeal cap |
| (1) Alexis O retractor |
| (3) 10mm sleeves |
| (1) 12mm sleeve |
| (1) Obturator for sleeves |
| (1) Instrument shield |
| (1) Introducer for the Alexis O retractor |



| Reorder No. | Alexis O Retractor Inner Ring Size | Qty/Box |
|-------------|------------------------------------|---------|
| C2A12 | 9.5cm | 1 |

Components

| |
|---|
| (1) GelSeal cap |
| (1) Alexis O retractor |
| (3) 10mm sleeves |
| (1) 12mm sleeve |
| (1) Obturator for sleeves |
| (1) Instrument shield |
| (1) Introducer for the Alexis O retractor |



| Reorder No. | Alexis O Retractor Inner Ring Size | Qty/Box |
|-------------|------------------------------------|---------|
| C2A15 | 7cm | 1 |

Components

| |
|---|
| (1) GelSeal cap |
| (1) Alexis O retractor |
| (3) 10mm sleeves |
| (1) 12mm sleeve |
| (1) Obturator for sleeves |
| (1) Instrument shield |
| (1) Introducer for the Alexis O retractor |

Visit [appliedmedical.com/gelpoint/v-path](https://www.appliedmedical.com/gelpoint/v-path) for more information.

The GelPOINT V-Path transvaginal access platform is used for gaining transvaginal access.

Please contact your Applied Medical representative for more information on availability. This information is intended for dissemination exclusively to healthcare professionals and is not intended to replace labeling and Instructions for Use (IFU). Please refer to the IFU for the indications, contraindications, warnings, precautions, instructions and other information.

© 2023 Applied Medical Resources Corporation, Rancho Santa Margarita, CA. All rights reserved. Applied Medical, the Applied Medical logo design, and marks designated with the symbol ® are trademarks of Applied Medical Resources Corporation, registered in one or more of the following countries: Australia, Canada, Japan, South Korea, the United States, the United Kingdom, and/or the European Union.